# COVADYE FLUORESCENT LABELLING KITS & SERVICES

**Custom Antibody Modification** 

From Research to Discovery

www.covalab.com



# **COVADYE** Simplify your labelling process

Labelling biomolecules are used to detect immediatly the presence and the location of your proteins or other biomolecules in biological systems (cells, tissus...).

Thanks to direct labelling, the amount of fluorescent signal detected is directly proportional to the amount of target protein in your sample. Consequently, no enzyme or binding compounds are required for the detection.

Covalab offers 2 solutions to increase your possibilities of coupling antibodies, proteins or peptides, with a wide range of fluorochromes :

### COVADYE LABELLING KITS

Taking advantage of a simple and fast labelling of antibody in your lab directly for particular wavelengths



## CUSTOM COVADYE FLUORESCENT LABELLING SERVICES

Benefit of a specific, precise and personalized labelling achieved by our technical expertise on a wide range of wavelengths





Antibody labelling kits offer an efficient and easy method for the fluorescent conjugation of your antibodies and molecules.

The kit provides everything required for the reaction and purification ensuring reproducibility and easy conjugate process.

In addition, to optimize your working time with a protocol ready to use, the kit is based on the principle of gel filtration to retain the labelled antibodies only.



#### **EFFICIENT**

RELIABLE

50 µg of antibody.

**READY TO USE** 

High coupling capacity decreasing the ratio of antibodies / fluorochromes.

Extremely stable conjugation able to label



#### **TIME SAVING**

Reproducible and optimized workflow achievable in 30 minutes before labelling incubation.



#### **EFFECTIVE**

High brightness and purity of the dye.



#### **COST SAVING**


One kit for all applications: Flow cytometry, Immunofluorescence, Western blotting...




### **APPLICATIONS**

#### **HIGH LABELLING CAPACITY**

Sample	DOL (Degree of Labelling)
Competitor 488	7, 7
CovaDye Micro 488 Labelling Kit	9,5

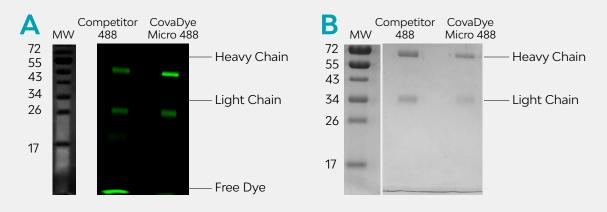


Electrophoresis Analysis comparing a competitor fluorescent dye and CovaDye Micro 488 Labelling Kit (kit5005). A : Revealed by Fluorescence 488 Filter B: Revealed by silver nitrate staining Sample loading (1 µg).

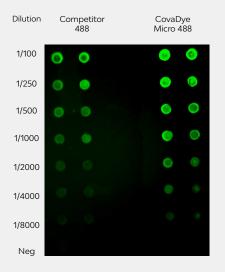

Dot Blotting analysis comparing a competitor fluorescent dye and CovaDye Micro 488 Labelling Kit (kit5005). Revealed by Fluorescence 488 filter. Mouse antibody spoted (1 µg) and anti-mouse (various dilutions).

#### STABLE AND SENSITIVE LABELLING

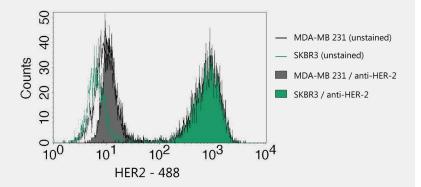

Flow cytometry analysis of cell surface HER-2 receptor using an anti-HER-2 antibody (1  $\mu$ g) directly labelled by CovaDye Micro 488 Labelling Kit (kit5005). Two cell lines were tested :

- MDA-MB 231 cell line, HER-2 negative (dash: unstained cells, grey: stained cells),

- SKBR3-cell line, HER-2 positive (dash: unstained cells, green: stained cells.



### **PRODUCTS**

#### TRIAL SIZE KIT FOR LABELLING 50-200 µG OF ANTIBODY

Perfect quantity for small amount of antibodies or for the optimization of your experimental conditions by testing different ratios antibodies / fluorochromes for saving your samples.

Products	Ex/Em (nm)	Packaging	ID Covalab	Code
CovaDye Micro 488 Labelling Kit	495-506	3 x 50 - 200 µg	kit5005	00118220
CovaDye Micro 550 Labelling Kit	550-565	3 x 50 - 200 µg	kit5006	00118222
CovaDye Micro 645 Labelling Kit	648-667	3 x 50 - 200 µg	kit5007	0018224
CovaDye Micro 770 Labelling Kit	774-790	3 x 50 - 200 µg	kit5008	00118226

#### PRODUCTION SIZE KIT FOR LABELLING UP TO 1 MG OF ANTIBODY

Products	Ex/Em (nm)	Packaging	ID Covalab	Code
CovaDye 488 Labelling Kit	495-506	1 x 1 mg	kit5001	00118209
		2 x 1 mg	kit5001	00118208
CovaDye 550 Labelling Kit	550-565	1 x 1 mg	kit5002	001182011
		2 x 1 mg	kit5002	001182012
CovaDye 645 Labelling Kit	648-667	1 x 1 mg	kit5003	001182014
		2 x 1 mg	kit5003	001182015
CovaDye 770 Labelling Kit	774-790	1 x 1 mg	kit5004	001182017
		2 x 1 mg	kit5004	001182018

# 2 CUSTOM COVADYE FLUORESCENT LABELLING SERVICES

#### LET US DO IT FOR YOU !

Direct labelling allows the use of several different fluorescent colors that help you in the study of major biological processes by giving you the ability to selectively identify your target.

Covalab offers the possibility of increasing your ability to label antibodies, proteins or peptides, with a wide range of fluorochromes.

#### **CUSTOM SERVICE**

Ultra-personalized services including technical support advice from immunology experts assisting you in your project.

#### INCREASE YOUR RESEARCH OPPORTUNITIES

High custom conjugation of a wide range of antibodies and other small molecules.

#### VARIETY OF FLUORESCENT DYE

Covering the broad range from **400 to 800 nm**.


#### **IDENTIFY ONLY YOUR TARGET**

Low non-specific backgound.

#### EFFECTIVE

One labelling for all applications: Flow cytometry, Immunofluorescence, Western blotting...

Label	Ex/Em	
CovaDye 425	425/476	
CovaDye 488	495/506	
CovaDye 530	529/561	
CovaDye 550	553/568	
CovaDye 570	573/612	
CovaDye 590	584/598	
CovaDye 630	629/646	
CovaDye 645	648/667	
CovaDye 678	677/703	
CovaDye 770	772/803	
CovaDye 800	802/825	
And many others		

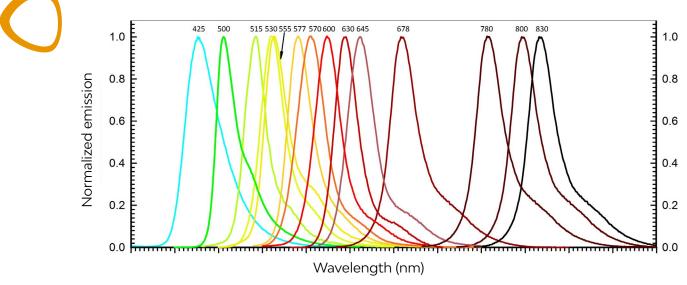

Emission Spectra of the CovaDye range

# SCIENTIFIC SUPPORT EXPERTISE HIGH QUALITY REACTIVITY - FXIBII **SINCE 1995**

(covalat

R&D in Biotechnology

#### **CONTACT US**

#### Covalab S.A.S.

11 avenue Albert Einstein 69100 Villeurbanne France

► +33 (0) 437 654 230
☑ enquiries@covalab.com

www.covalab.com